



Urea Broth Base (Diagnostic Stuarts Urea Broth Base)

M111

Intended Use:

Recommended for the identification of bacteria on the basis of urea utilization, specially for the differentiation of *Proteus*, *Salmonella* and *Shigella* species from clinical and non-clinical species.

Composition**

Ingredients	g/ L
Potassium dihydrogen phosphate	9.100
Di-potassium hydrogen phosphate	9.500
Yeast extract	0.100
Phenol red	0.010
Final pH (at 25°C)	6.8±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Suspend 18.71 grams in 950 ml purified/distilled water. Heat if necessary to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to 45-50°C. Aseptically add 50 ml of sterile U40 Supplement (FD048). Mix well and distribute in 10 ml amounts into sterile tubes.

Principle And Interpretation

Rustigian and Stuart developed Urea Broth (1). This medium is especially recommended for the differentiation of *Proteus* species from *Salmonella* and *Shigella* species in the enteric infection diagnosis (2), based on urea utilization (3,4). Gram-negative enteric bacilli are unable to utilize urea because of less nutrients and high buffering capacity of the medium. Urea Broth becomes alkaline as the utilization of urea by the organisms liberates ammonia during the incubation, indicated by pink red colour. All urea test media rely on the alkalinity formation and so they are not specific for urease testing. The utilization of proteins may raise the pH to alkalinity due to protein hydrolysis and excess of amino acids results in false-positive reaction. A medium without urea serves as negative control to rule out false positive results.

Type of specimen

Pure isolate

Specimen Collection and Handling

For food and dairy samples, follow appropriate techniques for sample collection and processing as per guidelines (5,6,7).

For water samples, follow appropriate techniques for sample collection, processing as per guidelines and local standards (8).

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (9,10).

After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions

In Vitro diagnostic use. For professional use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations

1. Prolonged incubation may cause alkaline reaction in the medium.
2. Also, all urea test media rely on the alkalinity formation and so they are not specific for determining the absolute rate of urease activity (4).

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Light yellow to light pink homogeneous free flowing powder

Colour and Clarity of prepared medium

Yellowish orange coloured clear solution in tubes.

Reaction

Reaction of basal medium (1.87gm in 95ml distilled water) at 25°C. pH : 6.8±0.2

pH

6.60-7.00

Cultural Response

Cultural characteristics observed on addition of sterile U40 Supplement (FD048) after an incubation at 35-37°C for 18-24 hours.

Organism	Inoculum (CFU)	Urease
<i>Escherichia coli</i> ATCC 25922 (00013*)	50-100	negative reaction, no change
# <i>Klebsiella aerogenes</i> ATCC 13048 (00175*)	50-100	negative reaction, no change
<i>Proteus mirabilis</i> ATCC 25933	50-100	positive reaction, cerise colour
\$ <i>Proteus hauseri</i> ATCC 13315	50-100	positive reaction, cerise colour
<i>Salmonella</i> Typhimurium ATCC 14028(00031*)	50-100	negative reaction, no change

Key : *Corresponding WDCM numbers.

Formerly known as *Enterobacter aerogenes*

\$ Formerly known as *Proteus vulgaris*

Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 2 - 8°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition Seal the container tightly after use. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (9,10).

Reference

1. Rustigian and Stuart, 1941, Proc. Soc. Exp. Biol. Med., 47:108.
2. Finegold and Baron, 1986, Bailey and Scotts Diagnostic Microbiology, 7th ed., The C.V. Mosby Co., St. Louis.
3. Christensen, 1946, J. Bact., 52:461.
4. MacFaddin J. F., 2000, Biochemical Tests for Identification of Medical Bacteria, 3rd Ed., Williams and Wilkins, Baltimore. Md.
5. American Public Health Association, Standard Methods for the Examination of Dairy Products, 1978, 14th Ed., Washington D.C.
6. Salfinger Y., and Tortorello M.L. Fifth (Ed.), 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
7. Wehr H. M. and Frank J. H., 2004, Standard Methods for the Microbiological Examination of Dairy Products, 17th Ed., APHA Inc., Washington, D.C.

- 8.Lipps WC, Braun-Howland EB, Baxter TE,eds. Standard methods for the Examination of Water and Wastewater, 24th ed. Washington DC:APHA Press; 2023.
- 9.Isenberg, H.D. Clinical Microbiology Procedures Handbook. 2nd Edition.
- 10.Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.

Revision : 00/2024

Packaging

M111-100G - Diagnostic Stuart's Urea Broth
Base (Urea Broth Base)


M111-500G - Diagnostic Stuart's Urea Broth
Base (Urea Broth Base)



HiMedia Laboratories Pvt. Limited,
Plot No.C-40, Road No.21Y, MIDC,
WagleIndustrial Area, Thane (W)
-400604, MS, India



In vitro diagnostic
medical device

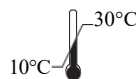

Storage temperature



CEpartner4U, Esdoornlaan 13,
3951DB Maarn, NL
www.cepartner4u.eu



CE Marking



Do not use if
package is damaged

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.